

FELINE HOUSING

R E I M A G I N E D



feline comfort

SUITE

INTRODUCING A NEW VISION FOR FELINE HOUSING

We are pleased to announce our new cat condos that are strong enough to withstand the harshest cleaning liquids and still have a silky, smooth surface for cat paws.



STRENGTH THROUGH TECHNOLOGY

We craft the condos from a revolutionary fine-celled foam board created with a new manufacturing process. This produces a hard but smooth outer structure that has a strong cellular core. Technically, it is an extruded PVC.

What that means to your feline clients is it is durable like outdoor signs and impervious to water and urine. Its strength can take a pounding from the largest cats without damage. The silky matte finish has just the slightest texture for secure footing.

We reached out to the construction industry to use the latest technology to create durable hinge that will stand up to repeated opening and closing. It is an injection-molded, fiberglass-reinforced nylon thermoplastic polymer.



LUXURY FELINE VETERINARY BOARDING

We've tried to get inside the feline brain through talks with our customers and thought leaders as well as time spent observing cats. We purposely reimagined how to make cats comfortable when spending time away from home in a kennel. The result is the **Feline Comfort Suite** that is easy on your eyes and soft on their paws.

Four door options are available:

The **clear, Polycarbonate door** covers the entrance while providing diagonal ventilation—a shape that makes cleaning easier. Two optional, adjustable vents in the back keep air flowing through the unit. All provides a clear view for cat and caretakers.

The **Tri-Panel** door combines a grill front with a polycarbonate band in the middle that gives the cat a bar-free view at eye level. The grill gives good ventilation through the unit.

Tri-Panel + Viewing Dome Door gives a special peeking area on the Tri-Door that increases feline interaction.

The **Full Grill Door** allows good airflow through the door and extra rods at the latch slide to protect paws.

Large litter space gives ample vertical room for a cat to take care of business in comfort. The litter has a **Grill Top with a Polycarbonate Bottom**. It comes with a clip for attaching information.

The litter features an optional **back vent**, similar to the two vents in the main unit. The door has **extra rods at the latch** to provide protection.

The **unique divider** has a cleaner edge that removes debris in the slot. You can "open" or "close" the portal status from outside the suite. **Optional storage** available.

Extra details include a **litter and fluid retention strip** at the litter front that helps avoid drips and a **lockable latch**.



U.S. Patent #7,805,808



LITTER AREA

(One Size)

Inside dimension:	16½"W x 28½"H x 28"D
Outside dimension:	18"W x 30"H x 28¾"D

LIVING AREA

(Two Sizes Available)

	Medium	Large
Inside dimension:	22½"W x 28½"H x 28"D	28½"W x 28½"H x 28"D
Outside dimension:	24"W x 30"H x 28¾"D	30"W x 30"H x 28¾"D

COMBINING OPTIONS FOR YOUR PERFECT FELINE COMFORT SUITE OR CONDO

You can create a custom-looking condo or suite assembly using our standard 24"W or 30"W main units and 18"W litter areas.

Note: The Tri-Panel + Viewing Dome can have the dome on the left or right side.

Main + Litter Door Styles



Tri-Panel



Tri-Panel + Viewing Dome



Polycarbonate



Full Grill



Standard Litter

Standard Condo + Suite Assemblies

These are a selection of the standard assembly units. They are also available with optional stationary leg or mobile platforms. Allow an additional 8½" of height for standard mobile platforms or 4" for height-reducing mobile platforms. Stationary platforms adjust to fit uneven floors, so plan for up to 9¼" of added height. **Note:** These assemblies all show pass throughs. Pass throughs are standard on Main + Litter units but optional on main unit (condo) only assemblies. These show Tri-Panel doors but all four door styles are available.



3½'W Cat Suite, Option A



4'W Cat Suite, Option B



5'W Cat Condo




7'W Cat Suite




7½'W Cat Condo



8'W Cat Condo

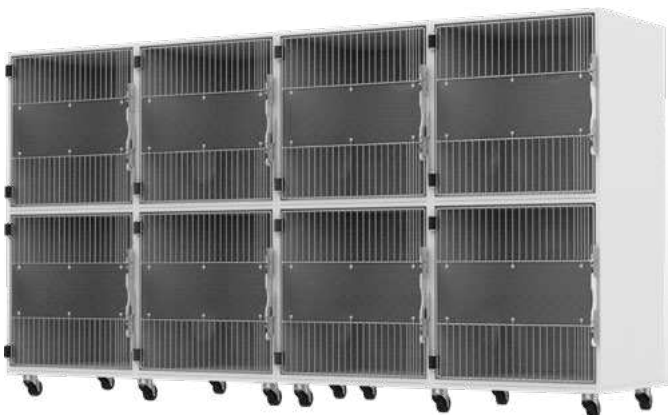

10'W Cat Condos

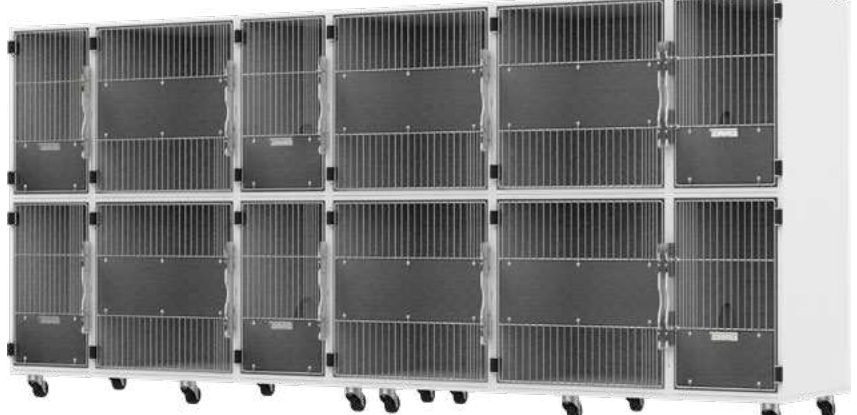

12'W Cat Suite

VENTILATION OPTIONS

Shor-Line offers two different ventilation options for Feline Comfort Suites and Condos, and both need to be specified at the time you place your order. Each requires holes to be cut in your unit.



PASSIVE VENTILATION

The main unit shown left features two 3¾" diameter vents placed high to encourage hot air to exit and circulate fresh air into the unit.



ACTIVE VENTILATION


The active ventilation system fits neatly and stylishly into the corner of the litter and/or main unit. It minimizes the reduction in floor space while allowing good air flow.

The longer unit creates a "tunnel" for sucking out the warmer air and venting it through a standard 4"-diameter opening at the top that attaches to your ventilation system. If you have a stacked assembly, the triangular unit fits snugly to the top corner and directs air up the longer tub and out your exhaust system.

Maintenance is simple with this system. Strong magnets hold the cover plates in place. Just flip them off and wipe clean. Then, replace it for a strong seal. A screwdriver allows you to adjust the vent opening and lock it down. Units are sold by the top and bottom piece.



Long Top Piece

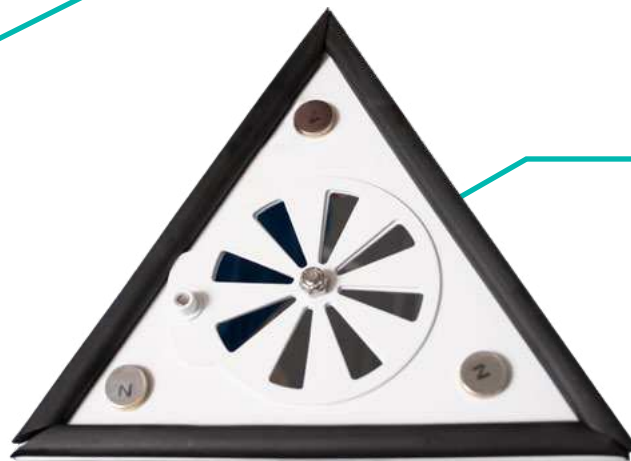

Triangular Bottom Piece

A NEW ERA OF FELINE BOARDING COMFORT HAS BEGUN

Our Feline Comfort Suite was developed with the latest feline housing and low stress design principles. Notice the many premium features that make the product appealing to cats as well as easier for caretakers.

We are pairing our craftsmanship with a revolutionary fine-celled foam board created with a new manufacturing process. It is technically referred to as an extruded PVC product. The result is a feline housing unit that could become as legendary and durable as our Stainless Steel cages. **Call TODAY to order an assembly!**



ACCESSORIZE

Optional Kat Kave: Two-level cave lets cats lounge on top, snuggle in the bottom or wind through the three openings below. Made of durable, easy-to-clean black, high molecular weight polyethylene.

SPECS: 20 $\frac{3}{4}$ "W x 14"H x 14 $\frac{3}{4}$ "D.

Privacy Quiet-Time Cover. Colorfast nylon cover snaps onto the cage to provide a space of privacy for the cat. It comes in blue and is machine washable. **Note:** Design for Feline Comfort Suites with Tri-Panel or grill fronts. Does not fit polycarbonate door.

SPECS: 12"W x 20 $\frac{1}{2}$ "H

